
♦lykk♦

KOMPRESSIOLAHKEET 2000

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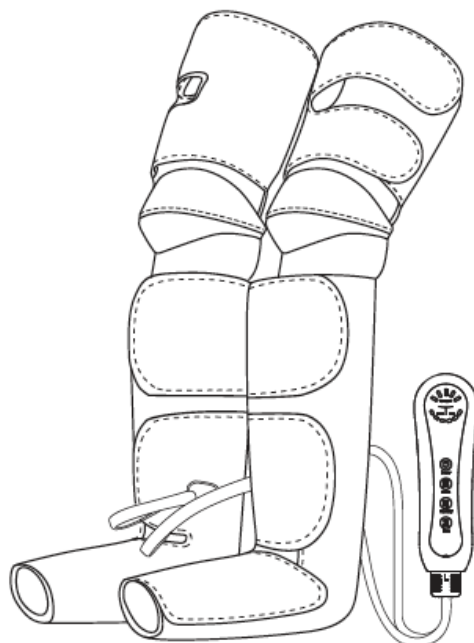
•lykke•

AIR COMPRESSION LEG MASSAGER 2000

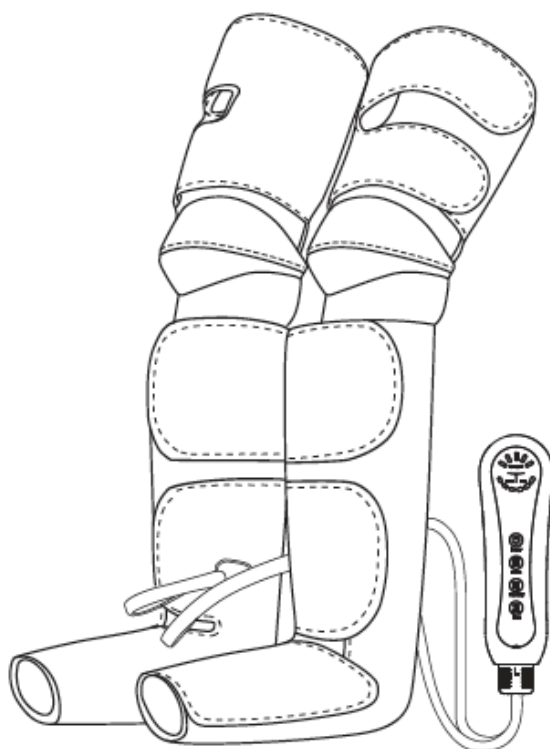IP-Agency Finland Oy - All rights reserved. Product and instructions for use are subject to change. Technical information is subject to change without notice.

IMPORTER:
IP-AGENCY FINLAND OY
HONKANUMMENTIE 13, 01380 VANTAA

Thank you for purchasing our products.

- Please read the instruction manual carefully before use, especially the "safety precautions", fully understand the scope of application and operation methods of the product, so that you can use it correctly and safely.

SAFETY PRECAUTIONS

Attention!: Do not use this device to supplement or replace any treatment you are currently receiving. If you encounter any unexplained pain, have received treatment for circulation or vascular problems, or use home medical equipment, please consult your doctor before use.

1. Before using our product, please be sure to fully understand the scope of application, operation methods, Precautions, etc.
2. If you have any product quality problems, please consult the seller. Do not disassemble modify or repair by yourself.
3. Please charge the controller before use.
4. This product is a household electronic health care product, which is not suitable for diagnosis, treatment and other commercial purposes, and cannot replace medicine and medical machinery.
5. Children and people with limited self-care ability and those with sensory defects are not allowed to use.
6. Use with caution for those with physical or mental disability or intellectual disability. Please use it under the Supervision of escorts. Any discomfort during use, please stop using it immediately.
7. Please use it properly in a safe environment.
8. If you feel unwell or have any discomfort during use, please stop using the product immediately.
9. Do not put this product in an unattended place, especially where children play around.
10. Do not use this product directly on areas with swollen or inflamed skin or open wounds.
11. Do not use in the bath or shower.
12. Do not place or store this product in a place where it may fall into the bathtub or sink. Do not immerse the product in water or other liquids.
13. Do not use pins or other metal fasteners on this device. Keep dry: Do not operate in a humid environment.
14. Do not operate in an environment with aerosol (spray) products or where concentrated oxygen is used.

Warning!

If you have any of the following condition please do not use this product.

- Pregnant or just after giving birth.
- People have a pacemaker or other medical implant.
- People suffering from heart disease or vascular disease.
- Abnormal blood pressure or receiving treatment to return blood pressure to a normal status.
- People suffering from fracture osteoporosis, spinal fracture fractures, muscle separation and other acute painful diseases.
- People suffering from cancer.
- People recently had surgery.
- People with sensitive skin or poor blood circulation.

If you encounter the situation below, please consult with your doctor before use



- People with malignant tumors (tumors).
- People with heart disorders.

- People who lose their sense of temperature.
- People with impaired consciousness.
- People with wounds at the treatment site.
- People who have a body temperature above 38 °C (with a hot period) (eg, acute inflammation Symptoms such as burnout, chills, changes in blood pressure, etc., and when they are weak) and who need to recuperate or when the physical condition is significantly worse.

PRODUCT INTRODUCTION

CONTROLLER

- 1: Pilot lamp
- 2: Power
- 3: Mode
- 4: Intensity
- 5: Heat
- 6: Leg Joint
- 7: Massage Air Tube
- 8: Leg Massager

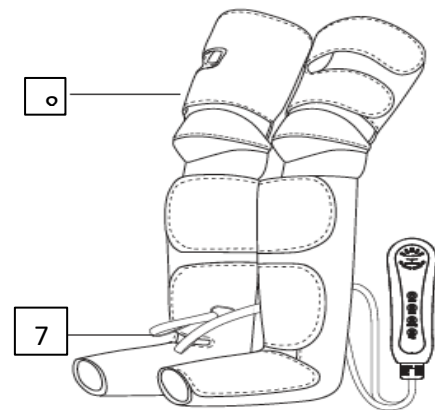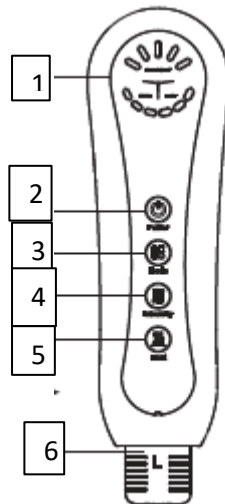PRODUCT SPECIFICATION

Product Name: Air Compression Leg Massager

Power Specifications: DC12V 2A

Power: 30W

Knee Temperature: 40/45/50 (Low-Medium-High level)

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Recommended using time: 15mins in each mode, 60 mins in one day

Massage Intensity: Foot: Weak: 30kpa Medium: 36kpa Strong: 38kpa

Calf: Weak: 22kpa Medium: 24kpa Strong: 30kpa

Thigh: Weak: 22kpa Medium: 24kpa Strong: 30kpa

Size Range of Each Part: Foot: about 26cm-36cm

Calf circumference: about 42cm-55cm

Thigh circumference: about 53cm-75.5cm

Product weight: About 1.65KG (3.4lbs)

Product Size: About 71.5x25x18cm (28.15" x 9.84" x 7.1")

USAGE METHODS


A. Check before use

-Check that the fabric is in good condition

- Check that it can work properly after powering on (please check this if it is not used for a long time).
- There may be signs of use on the skin. If you are worried about using marks, please wear on socks or pants before using the massager. It may leave marks on the calves feet, knees and thighs, and may not fade immediately, so please pay attention before going out.

NOTE: Do not use it on clothes and trousers with zippers or metal objects. It may damage the product.

B. How to connect the controller with massage hose?

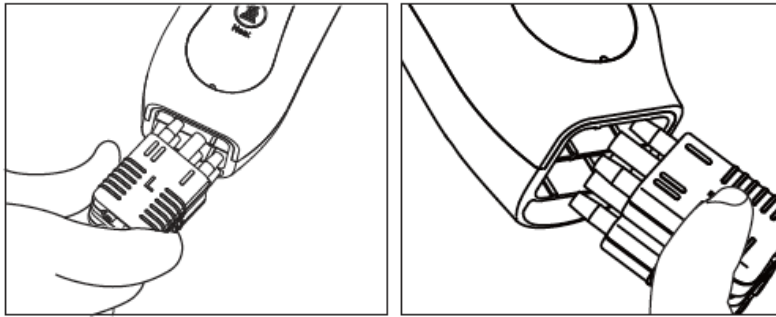
As shown above, insert the words "R" and "L" on the hose into the corresponding "R" and 'L' holes on the controller respectively.


Note: Make sure that the words on the hose are facing out, as we can see, that is, "L" is facing up, "R" is facing down. Otherwise the air pressure cannot be produced.

C. WEARING ORDER

1) Put the foot into the massager, start from the foot, fix the foot first, and then stick the velcro to the calf-thigh-knee in sequence, as shown below.
There are two velcros on the calf.

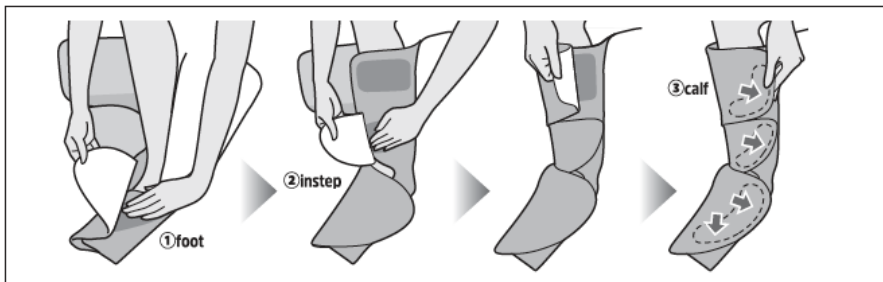- 2) There are 2 velcros on the thigh, which can be adjusted according to the individual's leg circumference.
- 3) After waiting for the installation of other velcro, it is easier to fix the velcro on the knee.

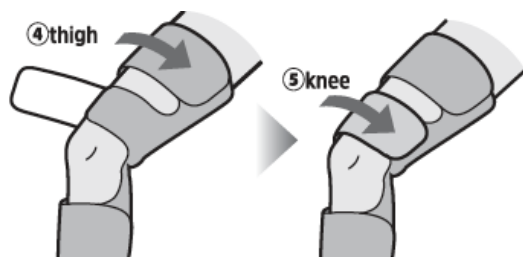Tips: It is easier to install if you sit on a chair. When installing, please be careful not to let Velcro stick to the clothes or stockings, especially sweaters, so as not to damage your clothes.


D. Ready to start using

POWER (Default timing is 15 minutes)

Short press power button to start the machine. The indicator light on.

Short press power button again to adjust the timing time. Each time you press it, it will increase by 5 minutes, up to 30 minutes. Long press for 3 seconds to stop the machine.

The screen displays:

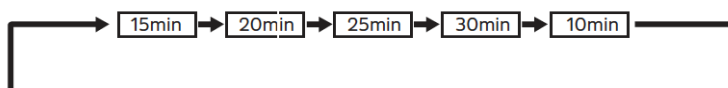Adjust the timing

[MODE] Press mode button to switch mode.

The sequence is as follows:

Sequence mode (Default power on mode M1)

In sequence mode, the upper compartment and the thigh are compressed and released while the lower compartment and the foot are compressed and released, and all four compartment air sacs are compressed and released simultaneously, there will be short rest period and then the cycle will start again.

Circulation Mode

In circulation mode, the upper calf chamber and thigh are compressed and maintained pressure is released. While the foot and lower calf chamber are compressed and maintained pressure is released, and the cycle procedure is repeated.

Combination Mode


In combination mode, combined sequence mode and circulation mode, the upper calf chamber and thigh were released after compression the foot and the lower calf chamber were released after compression, and all four chamber air sacs were released after compression, again circulating the procedure.

M1 Mode 1 Sequence

M2 Mode 2 Circulation

M3 Mode 3 Combination

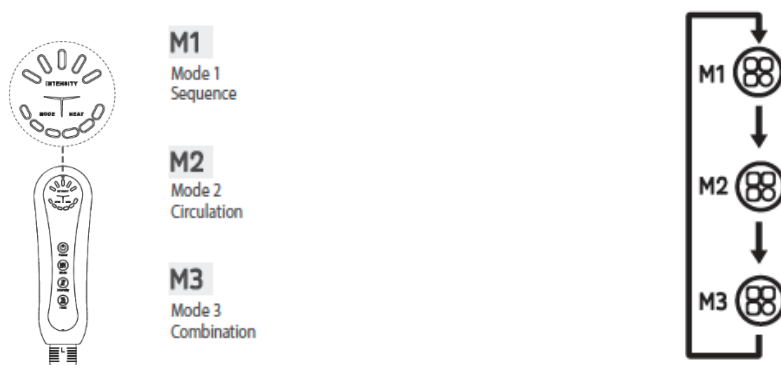(INTENSITY) Intensity Button

The sequence is as follows:

Lowest Intensity (Default power on Intensity L1)

Massage with low strength.

Second Lowest Intensity

Massage with second lowest strength.

Medium Intensity



Massage with medium strength.

Second Highest Intensity

Massage with second highest strength.

Highest Intensity

Massage with high strength.

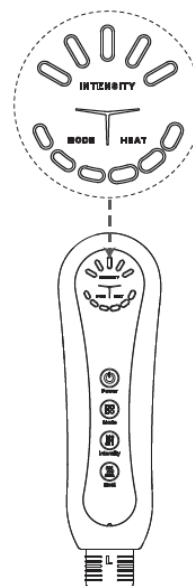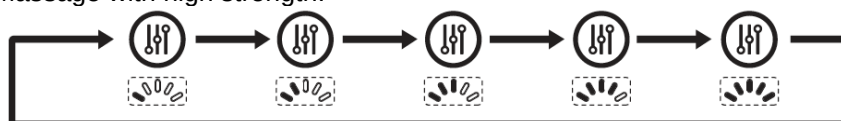(HEAT) Heater Button

The sequence is as follows:

Low Temperature (Default temperature is H1)

Massage with low temperature (about 40°C).

Medium Temperature


Massage with medium temperature (about 45°C).

High Temperature

Massage with high temperature (about 50°C).

Press heat button to switch the temperature mode,

When is displayed, the heating function is turned off.

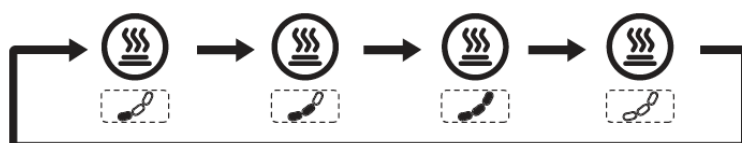NOTE: The product heats the kneet part.

Tips:

1. Fort he first use, it is recommended to start with the "weak" level of strength. You can increase the strength after getting used to it.
2. The air pressure is obvious after the massager is turned on, so when using it for the first time, install it according to the position of the velcro as much as possible. After getting used to It, adjust It according to your needs.
3. Please do not stand to use, sit down or lie down as much as possible to relax In order to better enjoy the massage.

CLEAN AND STORAGE

A.Clean

- * Make sure to turn off the power and unplug the hose.
- Never use thinner, benzine or alcohol (which may cause malfunction, broken parts, discoloration, etc.)

1) Wet the cloth and wring it out


(If you are worried about dust, please use a neutral detergent in the kitchen dilute to 5% with warm or clean water.)

2) Clean

Wipe the inside and outside of the leg massager remember not to wipe with large force.

* Note:


If rubbish or hair is caught in the hook and loop fastener, remove It with a toothpick, etc.
Do not Insert foreign objects Into the outlet of the hose plug.

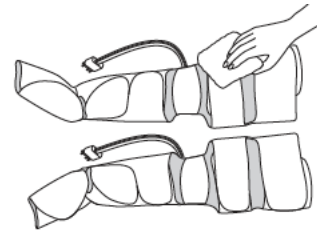B.Storage

1) Stick the Velcro on the feet, ankles and calves to the corresponding positions at once 1 2 3, and fold the foot part inward 4.

2) Fold the ankle and calf 5, round the connecting hose, and then fix it with a belt fastener 6 7.

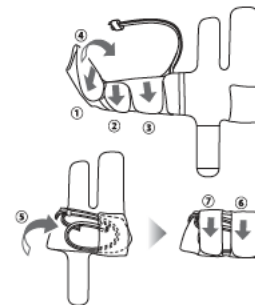C. Safekeeping

Put in a place where temperature and humidity do not rise and sunlight cannot reach.

1.Keep out of the reach of children.

2) Do not store other objects on the surface of the controller or massager (avoiding traces on the surface or color shift/ dlscoloratlon).

COMMON PROBLEMS

A. Under what circumstances should I not use it?

Answer:

- 1 During pregnancy or immediately after birth delivery.
- 2 Patients with severe blood circulation disorders cannot be used under the above circumstances. You can also see the details on page 2.

3 In addition, If you have a fever or any abnormality In your body, be sure to consult a doctor before use.

B. How many times can I use a day?

Answer:

For your safety,the massager will automatically shut down after 20 minutes of continuous use. We do not recommend that you use a massager for more than 1 hour. Prolonged use will over-stimulate muscles and nerves and may have adverse effects.

C. Can it be used for other parts than feet, calves, knees and thighs?

Answer:

It Is only for feet, calves, knees and thighs. Use In other parts may cause accident or Injury.

D. What size feet, calves and thighs are suitable for?

Answer:

Foot size: about 26cm to about 36cm.

Calf circumference: about 42cmt o about 55cm.

Thigh circumference: about 53cm to about 75.5cm.

E. The strength is too strong or too weak.

Answer:

The elasticity can be adjusted by the Velcro wrap method. The Velcro is loosely wound, the overall use strength and pressure will be relatively weakened. On the contrary,the Velcro can be adjusted properly, the overall intensity and pressure will be relatively strengthened.

** Velcro that stretches out after tightening may make the skin uncomfortable and can be covered with a thin handkerchief. Attach the handkerchief to the exposed Velcro.*

F. The air in the massager hose cannot exhaust.

Answer:

If a lot of air remains In the airbag, It may take some time to exhaust. Unplug the hose and wait for a while for the air to escape naturally.

G. The following sounds will be heard. Is it abnormal?

1. Velcro sound is crisp

This Is the sound of the Velcro slightly opening on the massager (this is the sound of the Velcro being pulled when the air enters the airbag).

2. Buzzing sound

This is the sound produced by the mode effect.

These are normal which is caused by the natural design of the massager structure.

TROUBLESHOOTING

Before requesting repairs, please confirm the following contents, and if errors persist after inspection, please contact your dealer.

Power cannot be turned on:

Check if the battery has power.

Stop working halfway:

1. Whether it is out of power.
2. Whether the hose plug Is loose.

3. Is the connecting hose broken?

Air does not enter

Is the connecting hose broken?

No temperature felt at the knee.

**Different feelings of warmth
between left and right knees**

1. The heating (heater) part will be heated in about 5 minutes.
2. The heated part does not fit closely to the knee, the gap is too large.
3. Whether the knee is in the set heating position.

Stress intensity is too weak:

1. Are the accessories loose? You can tightly wind the Velcro properly.
2. Do you wear thick socks or pants?

Stress intensity is too strong


Is the attachment too tight? You can relax the Velcro appropriately.

Air cannot be exhausted

Remove the accessory and pull out the hose plug to reinstall.

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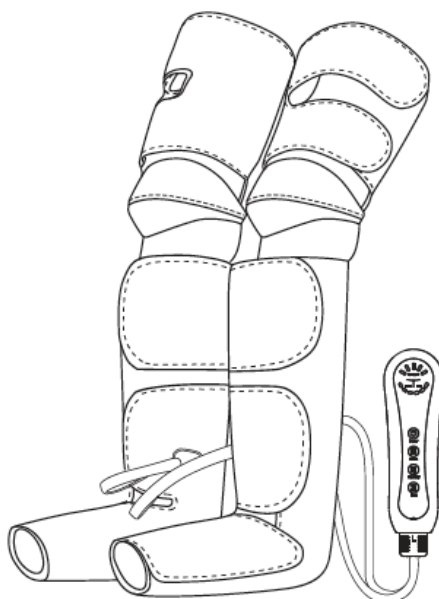Correct Disposal of this product	
	This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources.

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DECLARES THAT DESIGN AND MANUFACTURING OF THE FOLLOWING PRODUCT COMPLIES WITH THE FOLLOWING DIRECTIVES AND STANDARDS AND IS IDENTICAL TO THE PRODUCT.

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EMC Directive 2014/30/EU
RoHS Directive 2011/65/EU
LVD directive 2014/35/EU

Standards - Standardit – Standarder – Normen:

EN IEC 55014-1:2021, EN IEC 61000-3-2:2019+A1:2021; EN 61000-3-3:2013+A1:2019 +A2:2021; EN IEC 55014-2:2021, EN 60335-1:2012+A11:2014+A13:2017+A1:2019+A2:2019+A14:2019+A15:2021, EN IEC 60335-2-32:2021
EN 62233:2008, EN IEC 63000:2018

Place, date and signature:

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